





Over Wheelchair Table ASSEMBLY INSTRUCTIONS

- 1** Remove the contents of the box and check that you have all of the following pieces

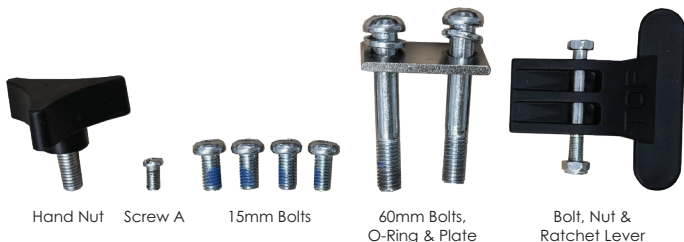


- 2** Fit the base underneath the mast, with the Split O-Ring/Washer between the plate and 60mm Bolt Head. Tighten with an Allen Wrench

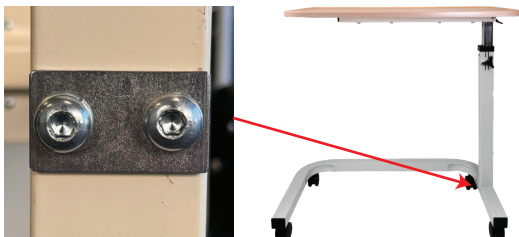


- 3** Attaching the table top to the mast:

- 3.1** Unscrew the bolt and nut from the ratchet lever
- 3.2** Slide the table into the mast
- 3.3** Attach table to mast using the four 15mm bolts
- 3.4** Re-assemble bolt, nut, and Ratchet Lever in mast

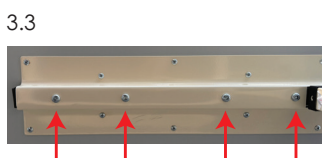

4

Remove the black locking screw in the mast and replace with Screw A

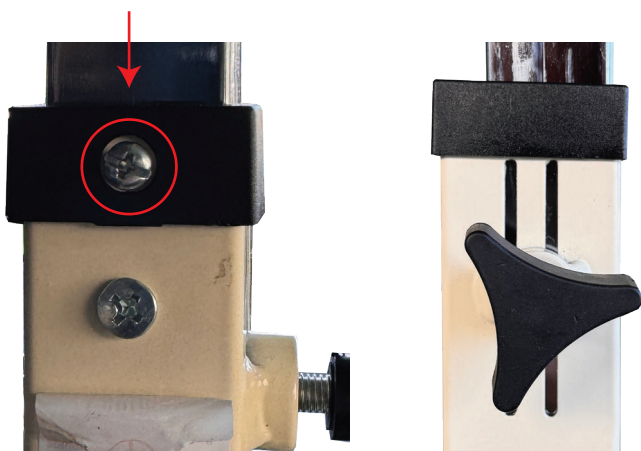

Note: When adjusting the height of the table, ALWAYS make sure that the Hand Nut is loose to avoid causing damage to the mast

5

See fully assembled table below. The table height can be adjusted by pulling the Ratchet Lever and raising/lowering the table top at the same time.



Note: Make sure to loosen the Hand Nut before adjusting, then tighten once at the preferred height to keep it in place.

CAUTION!

The Over Wheelchair Table is not designed to support heavy loads.

DO NOT sit on the table
DO NOT lean on the table for support when sitting down or standing up

Safe Working Load: 10kg

For more information visit www.goldfernmobility.co.nz